


### **3. sklop: ECMO - SISTEM ZA ODRASLE**

Skica za komplet cevi za levo ali desno podporo srcu s heparinsko oblogo

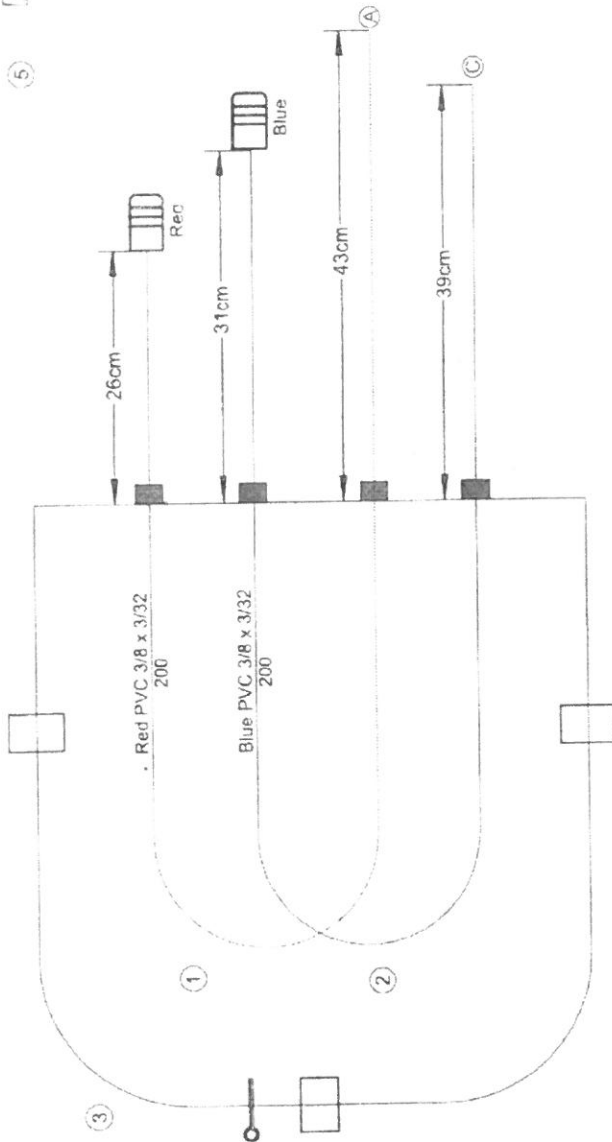

Breather Bag 860 x 450

Product Code / Revision Reference:

Description:

Page: 1  
Of: 1  
Date:

PACKAGING

Large Correx Pizza Box

Approved By:

A Remove connector on reservoir leaving 2cm of tubing  
Expand & bond 3/8" line into 1/2" reservoir inlet.

B Remove connector on reservoir leaving 3cm of tubing.  
Sleeve into 1/2" inlet a length of 3/8 x 3/32.  
Trim inlet to 2cm leaving a flush/flat finish.  
Internally bond prime line into inlet (no need for expanding).